

# APOLEX Tite

Exclusively Distributed by:



[www.technomed.net](http://www.technomed.net)

+91 85869 52121

[technomedindiapvtltd](#)



## Smart Lipolifting RF System

- ◀ Internal RF Tightening Probes
- ◀ RET Non-Invasive RF Tightening Probes
- ◀ Vaginal Tightening Probes



# The **Smart LipoLifting** World

LipoTite

SmartTite

MicroTite

SynergyTite

VaginalTite

## **Technologies**

Emusification

Lifting/Tightening

Smart LEDs for Temperature Indication

Energy Indications

Continuous/Pulse Mode

Energy Count

Real Time Temperature Indication

Portable and Light Weight

Thermo Sensor

# Lifting and natural skin stretching


Apolex Tite is a powerful **smart radio frequency device**. Its versatility and minimal invasiveness allow it to treat flaccidity in any part of the body, **emulsify fat** during liposuction, perform **facelifts** and **intimately rejuvenate** women.

Equipped with intradermal temperature sensors and an exclusive **smartLEDS** system.

Apolex Tite is positioned as one of the most revolutionary technologies of recent times.



Cannula with thermal sensors.

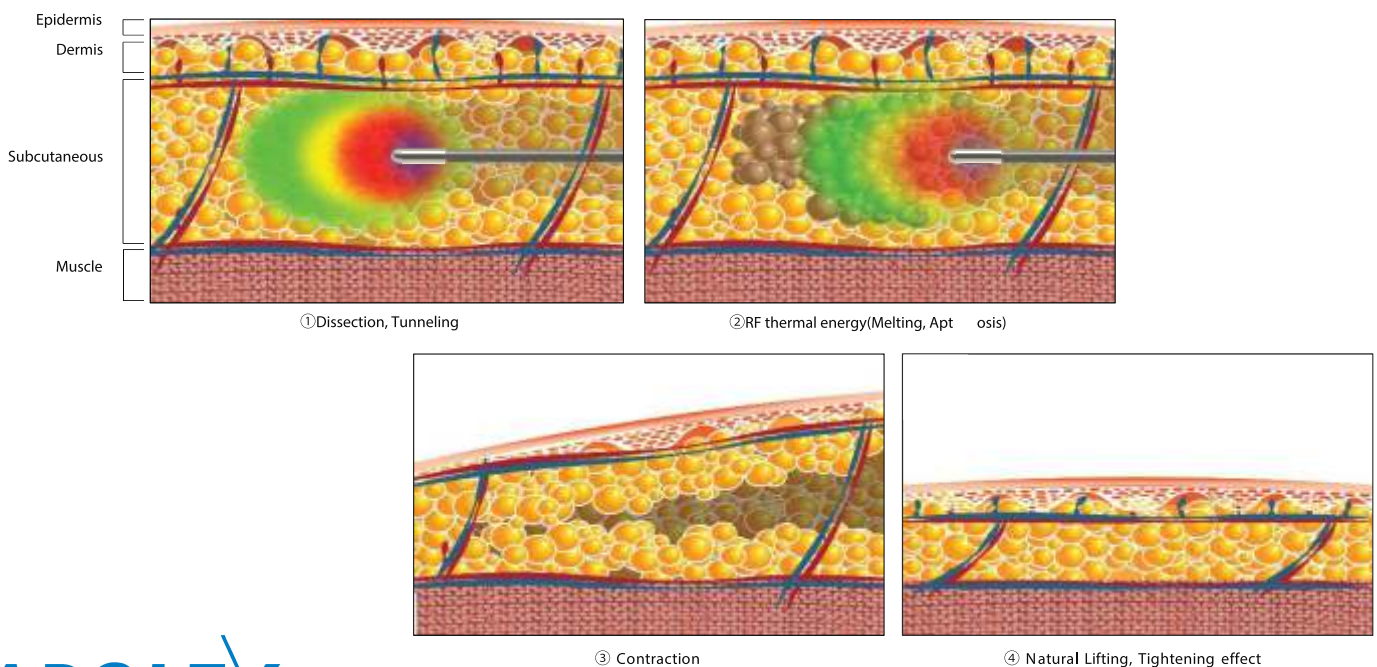

**APOLEX** Tite

# APOLEX Tite




Inset image details:

- 2.0Ø × 80mm
- 2.0Ø × 180mm
- ON/OFF button
- Cannula & Thermo Sensor
- 4.0Ø × 300mm
- LED indicator
  - Green : Ready
  - Blue : RF Operating
  - Red : Temperature Alert


Key features:

- Safety
- Effectiveness
- Economical efficiency
- Convenience

# APOLEX Tite Interface

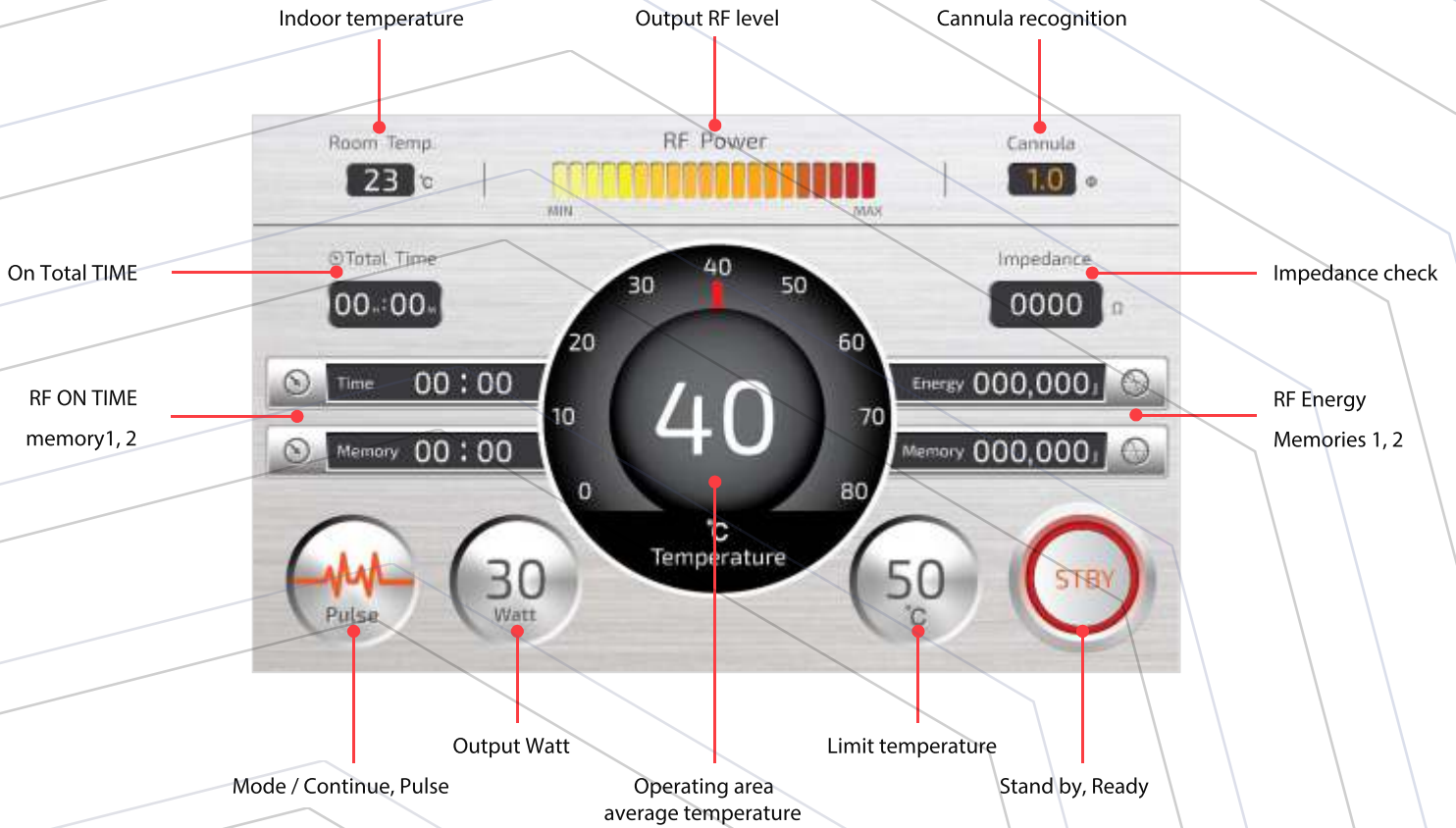

# APOLEX Tite Handpiece

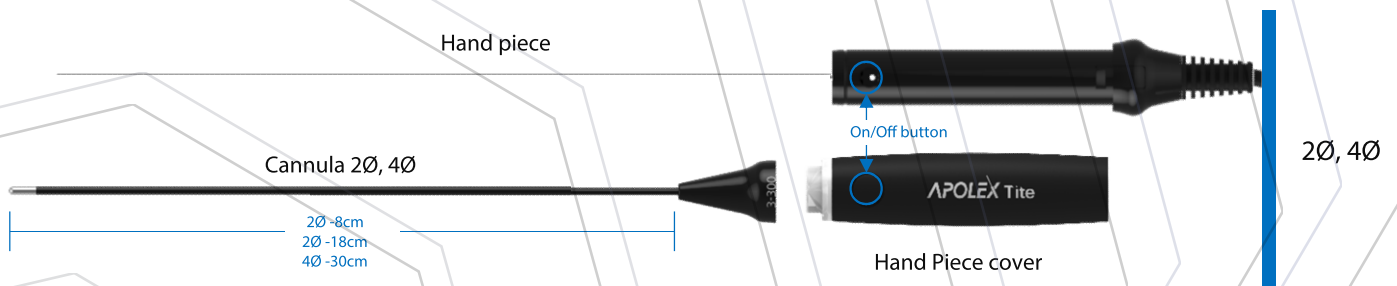

Hand piece  
Internal temperature, Impedance real-time control

Electrode Cannula  
Teflon insulated cannula

Hand piece Cover  
LED cover Lightning by temperature change (20, 40)

# LipoTite

The innovative **eMulsify Pulse** and its powerful **60W** energy emulsify adipose tissue for a more comfortable liposuction to the surgeon and safer for the patient causing internal heating contraction.


Its biological action to achieve retraction lies in a stimulus to the fibroblasts, it triggers new collagen and elastin fibers at the moment the energy comes into friction with the water molecules and increases the temperature for a prolonged time.

LipoTite and its technology pursues the idea of traditional plastic surgery cannulas, retaining their surgical ergonomics with a futuristic and highly intelligent design.

# LipoTite



300mm  
Ø4mm

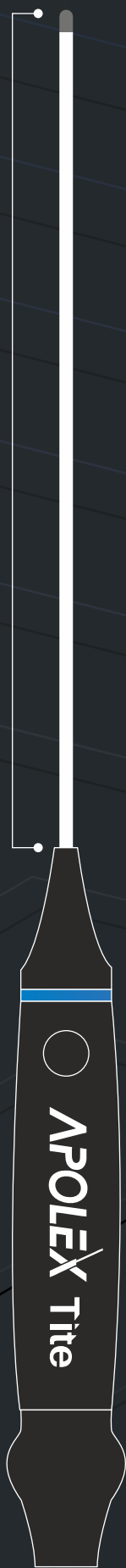

# SmartTite

SmartTite is a miracle of modern science, the first plastic surgery equipment in the world to develop last generation **smartLEDS technology**.

With a power of **60 Watts** and its **eMulsify** and **Tightening** Pulse technologies, **ApolexTite** emulsifies fat and retracts the skin in medium areas of the body such as arms, face (lower third, middle third) and neck.

SmartTite becomes a promise for the profiling of the mandibular area and jowls.



SmartTite



180mm  
Ø2mm



APOLEX Tite

# MicroTite

MicroTite is one of the smallest skin tightening handpieces in the world, measuring just **80mm** long.

Designed with ultra-light **sofTech technology**, allowing a comfortable approach without generating fatigue in the surgeon's hands.

Ideal for MicroTite small areas such as **nasogenian folds, eyelids** and **cellulite dimples**, it is also perfect for removing accumulated fat under the **armpits**.

# MicroTite



80mm  
Ø2mm

APOLEX Tite




# VaginalTite

VaginalTite is the most noble and revolutionary option for the **rejuvenation of the vaginal canal**, regenerating the **collagen** and **muscular elasticity** of the internal vaginal walls.

We can monitor the temperature during the treatment, offering safety to the operator and the patient.

Results with VaginalTite are predictable thanks to **EnergyCOUNT**; This technology counts the energy automatically, offering excellent benefits in sexual relations and urinary incontinence. The patients experience the changes from the first session and without any type of disability.

# VaginalTite

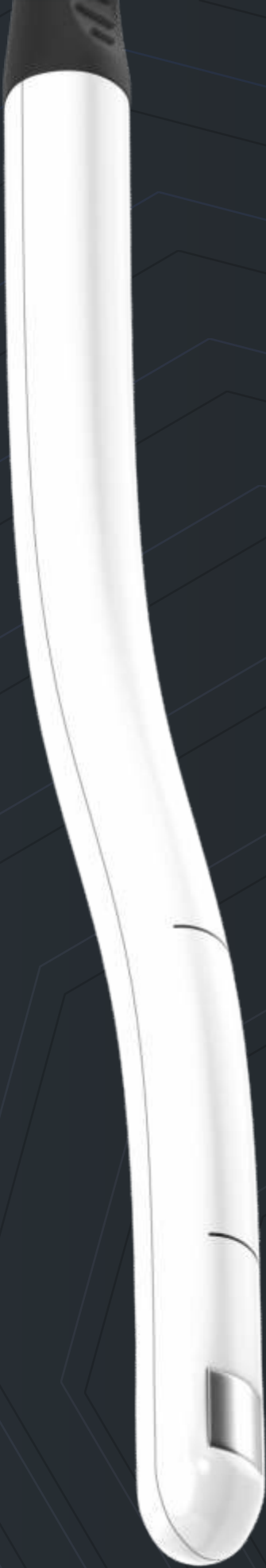

# SynergyTite

Betting on the optimization of abdominal marking results after **liposculpture**, **faceLIFT**, **bodyLIFT** and **fibrosis management**.


SynergyTite stands out as a leading technology in regeneration, **bodyshaping** and non-surgical procedures. Emitting energy deeply with a blast of electrons that come into contact with water molecules, triggering a process of reduction of fatty adipocytes.

It has been shown that this high energy delivery, at **460Khz**, improves blood oxygenation, favors microcirculation and weakens the connective walls that produce **cellulite**.


# SynergyTite




Ø30mm

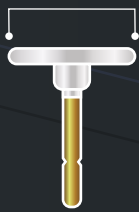

Ø50mm

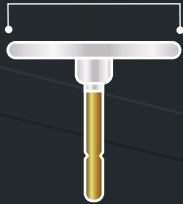

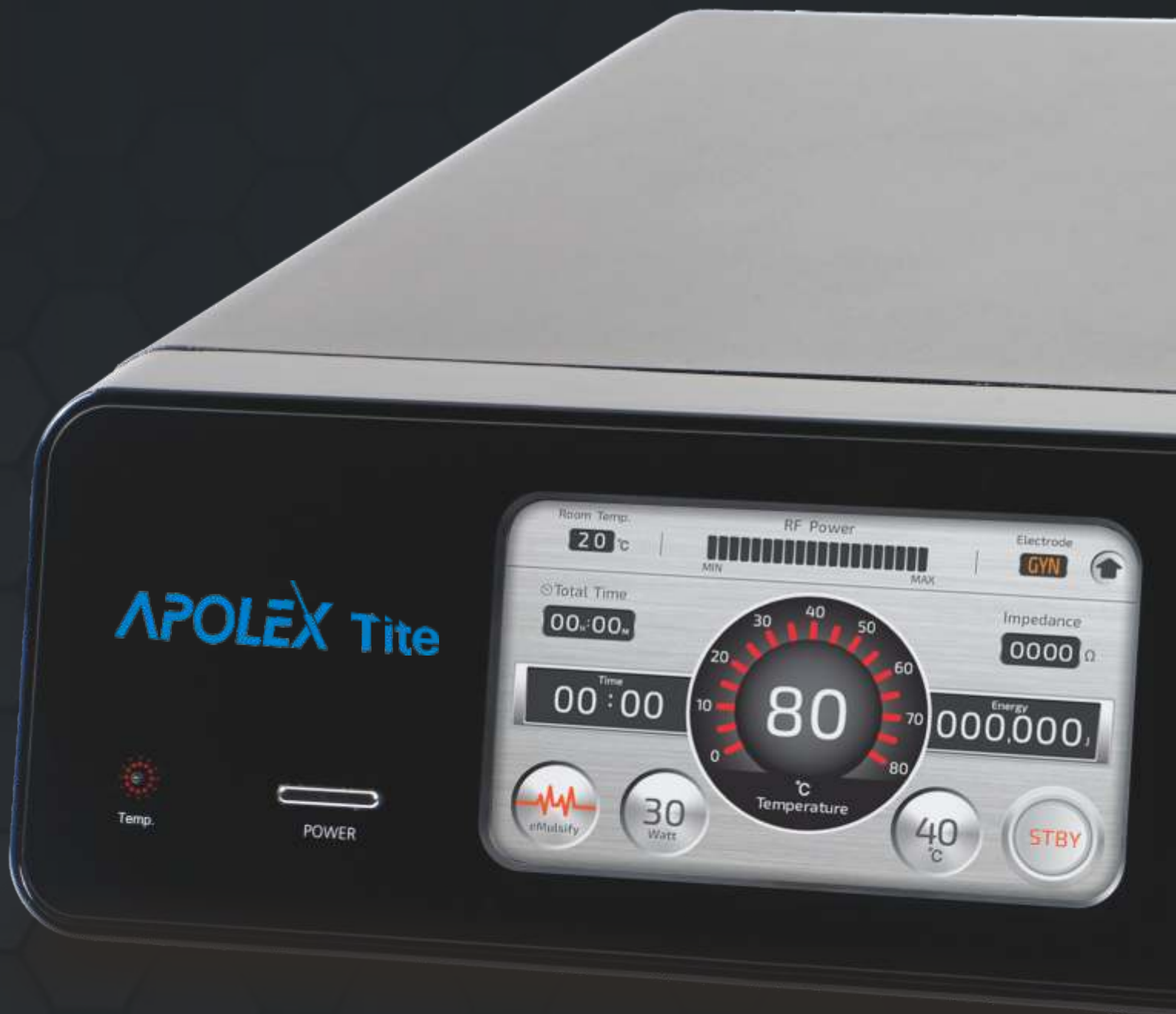

**APOLEX Tite**

Room Temp.  
**20 °C**



Electrode  
**GYN**

Total Time  
**00:00**

Impedance  
**0000 Ω**

Time  
**00:00**



Energy  
**000,000 J**

Temp.

POWER

**Multifly**

**30 Watt**

**40 °C**

**STBY**

**THE MODERN LIPO LIFTING**








# FEATURES

## smartLEDS

Patented technology that offers the operator operator **real-time** temperature monitoring, during the procedure in progress, without losing attention.

Patent No. 10-1380988.

State	Color	LED
Ready		Static
RF Operating		Static
Safe range		1s blink
Out of range		1/2s blink
Error		1/4s blink

## EnergyCOUNTER

Technology that counts the amount of energy delivered by each treatment zone. Offering predictable results and treatment protocols with high uniformity and precision.

## EnergyFLUX

FLX automatically measures the impedance from the conductive plate under the patient to the cannula.

## InnovativeDESIGN

Retraction has cannulas created for the modern surgeon, preserving its surgical ergonomics with a very smart and futuristic design.


APOLIX  
Tite

## Clinical cases

Retraction incorporates the most recent technology with excellent results:




**ABS**  
Courtesy: Dr. Humberto Uribe

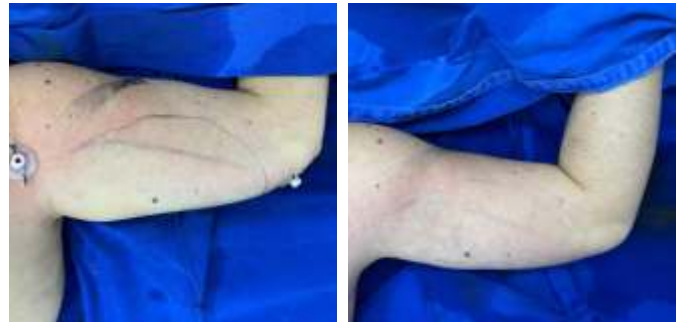

**ARMS**  
Courtesy: Dra. Sandra Ortega

## Clinical cases

Retraction incorporates the most recent technology with excellent results:

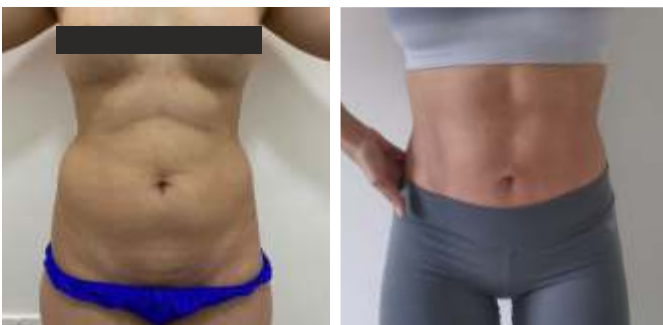

**ABS**  
Courtesy: Dr. Humberto Uribe



**ABS**  
Courtesy: Dr. Andrés Durán



## Datasheet

<b>Energy Supply</b>	AC 220V, 60Hz, 50VA
<b>Screen</b>	7" TFT LCD Touch Screen
<b>Dimensions</b>	392 (wide) x 414 (depth) x 123 (high) mm
<b>Weight</b>	7.22kg / 15.9lb
<b>Output Power</b>	50w
<b>Output RF</b>	1~50w
<b>Pulse</b>	Continuous / Pulse
<b>Output Frequency</b>	480KHz
<b>Handpieces</b>	3.0Pi / 2.0Pi / 0.8Pi / VA

# APOLEX Tite



## SMART LIPOLIFTING with Liposuction

### RELATED PRODUCTS



**MICROAIR<sup>®</sup> PAL**  
POWER ASSISTED LIPOSCULPTOR  
*For Surgery. For Life.™*



**LipoSaver**  
Ultrasound Assisted Liposculptor



**MEDCO ASPIRATOR IV**



**FAT COLLECTOR**



**TECH-3000 INFILTRATION PUMP**

Manufactured by:

**CHUNGWOO**

**CHUNGWOO MEDICAL**

(Gasan-dong, Woolim Lionsvalley-2), 614, 2,  
Gasandigital1-RO, Geumcheon-gu, Seoul, Korea  
Call us +82-2-2027-2200  
Mail : marketing@mycw.co.kr

Imported & Marketed by:

**Technomed India Pvt. Ltd.**

Regd. Office : C-19, DSII DC Complex,  
Kirti Nagar, New Delhi-110015

+91 8586952121 | Office : 011 - 41420555

E-mail : info@technomed.net | technomedindia@outlook.com

Website : www.technomed.net | www.pressuregarment.com