
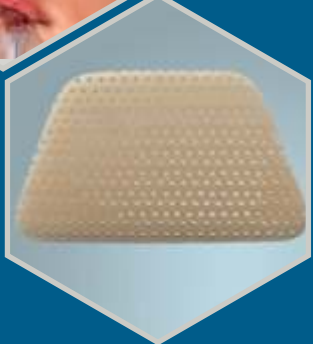
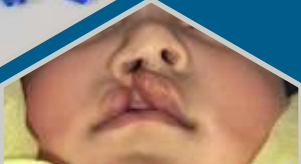



NOSTRIL RETAINER, INTERNAL & EXTERNAL NASAL SPLINTS

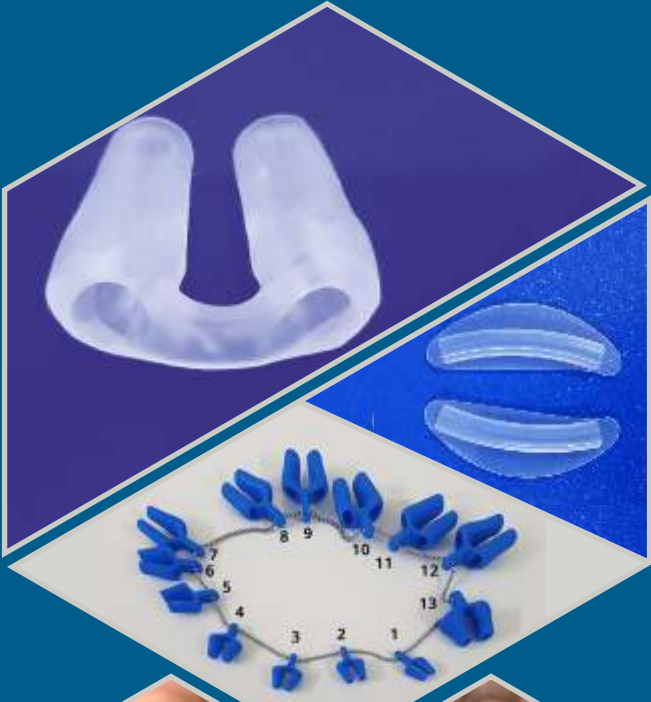

www.technomed.net
+91 85869 52121

Features

- This product prevents the deformation of the nasal foramen after the surgery for cleft lip or palate, rhinoplasty, or augmentation rhinoplasty without imposing a burden on the patient.
- It is available in 13 sizes each for the standard and long types so that it can be used in a size appropriate to each patient.
- It is made of soft silicone.

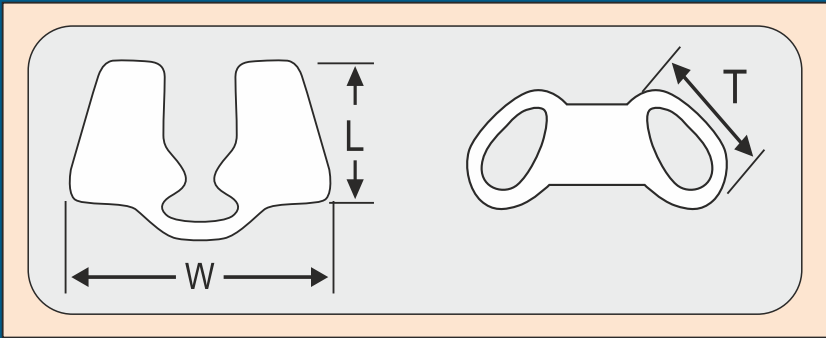

No	Device Dimensions W*L*t (mm)	Device Model
1	13*14*6	Pediatric
2	15*16*7	Pediatric
3	17*17*8	Pediatric
4	19*18*9	Pediatric
5	21*18*10	Pediatric
6	22*20*11	Pediatric
7	22*23*11	Adult
8	24*23*12	Adult
9	26*23*13	Adult
10	28*23*14	Adult
11	29*24*15	Adult
12	31*24*16	Adult
13	32*24*17	Adult
14	33*24*18	Adult

Internal Nasal Splints for Septoplasty

- The Internal Airway Nasal Splint, commonly known as the "Doyle Splint", provides support to the septum while maintaining an open airway for natural air flow after septal surgery.
- Internal splints are usually used after a septoplasty, or correction of a deviated septum. They're made up of two pieces of plastic or other firm but flexible material inserted inside the nose to support the bridge between the nostrils, or septum
- This intranasal splint reduces the risk of the airway closing from accumulation of blood and mucus by allowing the lumen to remain open. The Internal Airway splint is made of soft medical grade silicone to aid with insertion and removal and includes suture holes for suture placement.

- Sterile, single-use
- Increase airflow
- Soft medical-grade silicone





Surgiform[®]
FACIAL IMPLANTS



PIEZOTOME M+
Piezo Ultrasound Surgery Unit
SOLO



www.technomed.net

+91 85869 52121

TECHNOMED INDIA PVT. LTD.

Regd. Office : H-107, Shivaji Park, Punjabi Bagh, New Delhi-110026

+91 8586952121 | Office : 011 - 41420555

E-mail : info@technomed.net | technomedindia@outlook.com

Website : www.technomed.net | www.pressuregarment.com